






Independent 4 Life

£325.55 Inc. VAT

# AKW Advanced Wall Mounted Extra Wide Bariatric Fold-up Moulded Blue Padded Seat w/ Support Legs

Product Code: 04110P

## Product Images

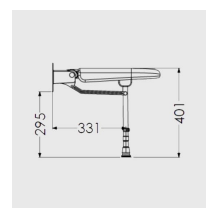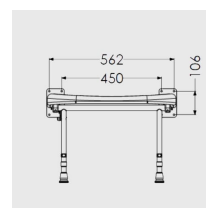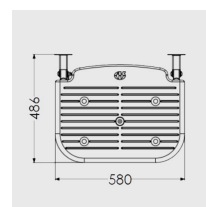



## Description

---

AKW 4000 series extra-wide bariatric fold-up padded shower seat with support legs, white moulded with inset blue padding, designed for accommodating varied individual needs. Extra wide seat and padding provide additional comfort for bariatric users, contrasting blue coloured padding and white mould aid the visually impaired. The seat pad can be easily removed making cleaning easier.

At the point of installation, seat height can be adjusted between 390mm – 640mm to maximise comfort by accommodating the clients' needs.

Can be folded away and stored against the wall when not in use, allowing the shower to be accessible to all.

This AKW shower seat has a maximum load capacity of 254kg / 40 Stone, wall mounted and is designed for easy cleaning.

### Key benefits:

- Made from stainless steel, ensuring durability for long-term use
- Comfortable in use
- +/- 13mm fine adjustment to support leg, allowing for gradient in floor to optimise comfort on a level seating
- Seat pad is easily removable for cleaning
- Folds compactly against the wall when not in use
- Ideal for any situation where side transfer is required
- Cushion fixing lugs prevent the cushion from moving during side transfer

### Technical details:

- Maximum load capacity (kg / Stone): 254 / 40
- Projection when folded (mm): 155
- Projection of seat from wall (unfolded) (mm): 486
- Width of seat (mm): 580
- Maximum height (mm): 640

- Minimum height (mm): 390
- Materials: Rust-proof stainless steel

**Warranty:**

- Lifetime manufacturer warranty

**Note:**

- The lifetime warranty excludes normal wear and tear in general domestic applications and in accordance with AKW installation instructions. It applies only to products being installed for the first time and according to relevant local and national regulations
- It is the responsibility of the installer to install the seat in accordance with best practice
- Support legs must not exceed a maximum of 1 in 40 shower base slopes when placed in the down position. Curved surfaces such as corner of bathtubs are not acceptable
- Wall fixings are not included

## Additional Information

---

Manufacturer	AKW / Contour
Operation	Manual Folding
Configuration	Seat Only
Seat Width (mm)	500 - 599
Height Adjustment (mm)	0 (Fixed)
Folding Seat	Yes
Padded	Yes